

A090英文说明书

单页尺寸：120mm*85mm



Far Infrared Vibration Heat Patch
HXR-A090 User's Manual

Danger

Before using this device, please carefully read the manual and comply with the operating guidelines of the "Instructional Manual", the company will not responsible for the consequences of improper use or non-compliance with the given instructions for the correct use of particular product.

This product is a personal care massage device and cannot be used as a substitute for medicines and medical device.

Caution

- Do not use this device near to the water or moist nor to the flammable.
- Do not use this device in the high electromagnetic environment, high temperature environment or in humid conditions.
- Do not use this device in the flammable and explosive places e.g gas stations, high dust workshops etc.
- Do not disassemble the device by yourself. if the device failure occurs, please contact the after-sales service to deal with the failure or to do its maintenance.
- Do not plug-in or un-plug the USB charging cable with the wet hands in order to avoid the danger of electric shock or any other kind of damage.
- Do not use this device with the other medical devices. Do not use more than 1 units of this product at the same time, otherwise it may cause the physical discomfort.
- Do not use it while driving a car, as it may cause physical discomfort due to excessive stimulation and that might lead to an accident.
- Do not use for any other purpose other than the massage as it is clearly stated in the instruction manual to avoid accidents and any kind of malfunctioning.
- If any kind of discomfort is caused by the excessive stimulation during its use, Please stop using the device immediately.
- This device should not be used by people those who are physically incompetent, have sensory or mental impairment, lack of relevant experience.
- Please do not give this device to the children or people without the ability to express consciousness to avoid accidents or physical discomfort.
- This device has heat, vibration and infrared red light function, please use with caution if you are not sensitive to the high temperature, vibration and Infrared light intensity.
- Do not use the device in the following conditions: pregnant women, high fever, malignant tumors, untreated heart surgery, acute illnesses, or other physical diseases.
- This device contains lithium battery, please do not disassemble, do not modify and do not throw into fire in order to avoid danger.
- When disposing of this device, please dispose off it in accordance with the relevant environmental protection regulations of the city.

Product parameters

Device Name	Far Infrared Vibration Heat Patch	Device Model	HXR-A090
Device Size	800mm*110mm*35mm	Rated Power	18W
Battery Capacity	3.7V 1900mAh	Charging Voltage	5V, 2A
Device Weight	161g	Charging Port	USB Type-C

*Note: Because the belt has good elasticity, so there may be ±5CM error between the measurement size and the actual size.

What is included in the package ?



Main Device*1



Type-C Charging Cable*1




Instructions Manual*1




Bundling Belt*1

Functions and names of each part

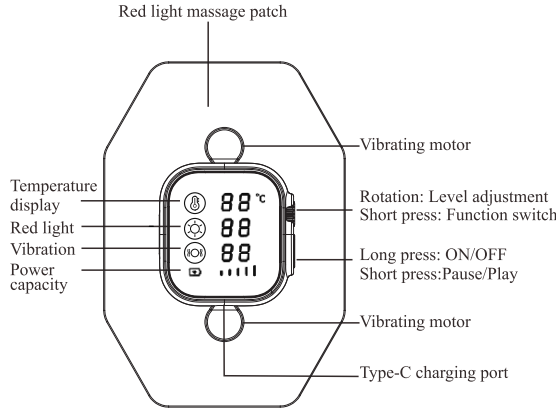

Device operating instructions


Control box panel

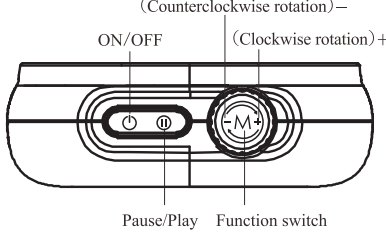

- Long press button to turn on, the LED light below the button will be turned on into the standby mode, if there is no other operation within 60 seconds, the device will be automatically turned off.
- During the Standby mode, then control the round knob to operate the device, the default working time is 15 minutes, so after 15 minutes the device will be automatically turned off, if the user need to continue to use, then the user need to restart the device.
- **Round knob** (Turn it clockwise to increase the gear. Turn it counterclockwise to decrease the gear.) :
 1. Press once temperature icon flashes, then rotate to select temperature, the adjustable range of the hot compress gear is: 38°C-48°C. (Actual temperatures may vary by ±3°C due to environmental factors, please adjust it according to your actual body feeling).
 2. Press again red light icon flashes, then rotate to select the light intensity, infrared light totally has 3 levels.
 3. Press again vibration icon flashes, rotate to select vibration intensity, vibration mode totally has 3 levels.
- **Pause button**: Short press once during the working process to suspend the function, and then short press to start the function.
- Long press button to turn off the device.

Assembly and wearing method

1. Pass the host through the center of the strap, as shown in figure ① below
2. After the installation is complete, as shown in figure ② below
3. After assembly, wrap the product around the part to be used and secure it with the velcro at the end of the strap to start using.

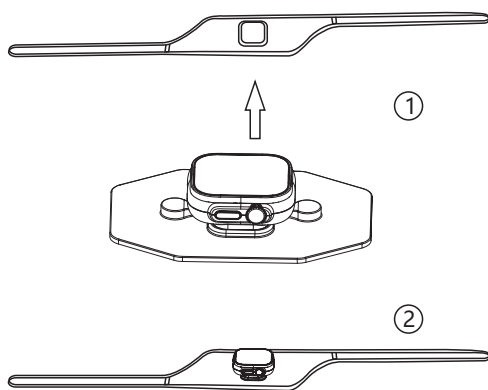

Applicable body parts



Neck



Arm/Leg



Wrist



Knee/Elbow




Ankle

Device charging

If it is used for the first time or has not been used for a long time, please charge it as follows.

1. Plug the Type-C port of the USB cable into the charging jack of the product.
2. Insert the USB cable into the power adapter and plug-in the power adaptor into the power socket.

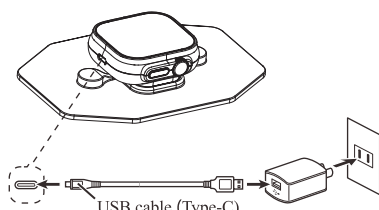

3. After charging starts, the power indicator of the host LED display blinks
4. The product takes about 2 hours to fully charge. After it is fully charged, the powerindicator light will be on for 5 bars. After charging for another five minutes, the atmosphere light will go out. Please unplug the power cord in time after it is fully charged.

*Note: Do not use this device while charging.

Troubleshooting

Possible conditions	Possible causes/handling methods
Cannot start after pressing the power button	Please check whether the device has the sufficient power or not Please confirm whether the device start-up operation is standardized (press and hold the on/off button for more than 1.5 seconds)
Stopped during its use	The device will automatically stop working after 15min of use, which is a normal phenomenon, if you need to continue to use the device then need to turn on the device again. Please confirm whether the device is in low battery status, if the battery power is insufficient, please charge it and use it again after its fully charged.
Unable to charge or very short life span	If the ambient temperature is particularly high or low, the device will not start, please charge at 0-45°C. Please check whether the USB connection cable and the host power adapter are in good contact The first charge or long time no use, please first full charge before its use. The device cannot be charged during startup, please turn off the power before charging

*Note: If the problem still exists after following the above steps, please cut off the power and contact customer service to deal with it.

Maintenance and storage methods

Maintenance method

1. Be sure to unplug the power cable before its maintenance.
2. Do not use detergent, thinner, gasoline or other chemicals to wipe or clean the device.
3. When the device is dirty, it can be wiped with a dry soft cloth or dry towel to avoid washing in running water or other liquids.
4. Please do not disassemble, repair or modify this device without the permission to avoid accidents or any kind of other possible failures.

Storage method

1. Be sure to unplug the power supply during storage.
2. Please do not keep it in a damp place or a place exposed to direct sunlight.
3. If the device is not used for a long time, please fully charge the device and put it in the storage bag, then put it in a dry and ventilated place.
4. Do not store it in a place where infants or pets can touch it.
5. Please do not place it on cloth, paper products or wooden furniture that are easy to fade.