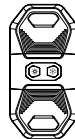
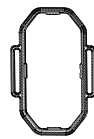
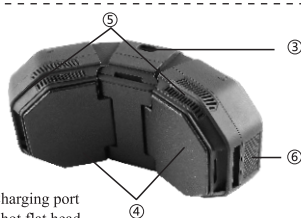
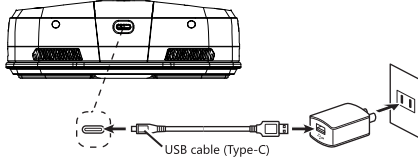


HXR-A130英文说明书八折页双面铜版纸彩色印刷

单页尺寸：70*90mm

HXR-A130  Recovery Therm Cube Instructional Manual	⚠ Caution ● Do not use this device near to the water or near to the place where there is fire ● Do not use the device in high electromagnetic environment, high temperature environment or humid environment. ● Do not use the device in flammable and explosive places such as gas stations, high dust workshops and gas stations. ● Do not drop or disassemble the device by yourself. Please contact the after-sales service to deal with the situation if the device failure occurs. ● Do not plug and unplug the USB charging cable with wet hands to avoid electric shock or accidents. ● Do not use this device with other medical devices or more than 2 units of this device, otherwise it may cause physical discomfort. ● Do not use it while driving a car, as it may cause physical discomfort due to excessive stimulation and that might lead to an accident. ● Do not use for any other purposes other than massage as stated in the instruction manual to avoid accidents and malfunction. ● If discomfort is caused by excessive stimulation during use, stop using the device immediately.	● This product is intended for general massage and physical comfort purposes only, not a medical device, and cannot replace professional medical care. ● This device should not be used by people those who are physically incompetent, have sensory or mental impairment, or lack relevant experience and knowledge. ● Please do not give this device to the children or people without the ability to express consciousness to avoid accidents or physical discomfort. ● This device has heating and cooling functions, so please use with caution for people who are sensitive to temperature change or who are allergic to heat or cold. ● Do not use the device in the following cases: pregnant women or people with high fever, malignant tumors, untreated heart surgery, acute illnesses, or other diseases. ● This device contains lithium batteries, please do not disassemble, modify or throw into the fire to avoid danger. ● When disposing of this device, please dispose of it in accordance with the relevant environmental protection regulations of the city, otherwise it may cause environmental pollution.	Product parameters <table><tr><td>Product Name</td><td>Recovery Therm Cube</td><td>Product Model</td><td>HXR-A130</td></tr><tr><td>Product size</td><td>149mm*76mm*50mm</td><td>Rated Power</td><td>20W</td></tr><tr><td>Battery capacity</td><td>3200mAh</td><td>Charging Voltage</td><td>5V</td></tr><tr><td>Product Weight</td><td>387g</td><td>Charging Port</td><td>USB Type-C</td></tr></table> Product list  Main Device*1  Strap Frame*1  Type-C Charging Cable* 1  Instructions*1  Strap*3	Product Name	Recovery Therm Cube	Product Model	HXR-A130	Product size	149mm*76mm*50mm	Rated Power	20W	Battery capacity	3200mAh	Charging Voltage	5V	Product Weight	387g	Charging Port	USB Type-C	Function and name of each part  ① Cool function Button ② Heat function Button  ③ Type-C Charging port ④ Cold and hot flat head ⑤ Air inlet (working in cooling mode) ⑥ Air outlet (working in cooling mode) Note: When the device is in cooling mode, please make sure that ⑤ air inlet ⑥ air outlet is not blocked to avoid affecting the cooling effect of the device.	Product operating instructions Turning On and Adjusting the Heat Function ● When the device is off, long - press the Heat Function Button for 3 seconds to turn it on. The device will default to Heat Function Level 1, and the red indicator light will illuminate. ● There are 3 levels for the heat function, with temperatures: 38±3°C, 40±3°C, and 42±3°C. Short - press the Heat Function Button to switch levels. The number of red indicator lights corresponds to the current level (1 light for Level 1, 2 lights for Level 2, 3 lights for Level 3). Turning On and Adjusting the Cool Function ● When the device is off, long - press the Cool Function Button for 3 seconds to turn it on. The device will default to Cool Function Level 1, and the blue indicator light will illuminate. ● There are 3 levels for the cool function, with temperatures: 15±3°C, 10±3°C, and 5±3°C. Short - press the Cool Function Button to switch levels. The number of blue indicator lights corresponds to the current level (1 light for Level 1, 2 lights for Level 2, 3 lights for Level 3).	Turning On and Adjusting the Cool & Heat Cycle Mode ● When the device is off, press and hold both the Cooling and Heating buttons for 3 seconds to enter the automatic cycle mode, and the indicator light will flash blue. ● In this mode, the cool and heat modes switch automatically every 3 minutes, starting from Cool Function Level 2 by default. During the cool or heat phase, you can short - press the corresponding (cool/heat) function button to adjust the level under the current mode. Shutdown Logic ● In all modes, the device will automatically shut down after continuous operation for 15 minutes. During operation, long - press any function button for 1 second to manually turn off the device. If you need to continue using it, you must turn on the device again. *Note: The above cool temperatures are tested at an ambient temperature of 20°C. During actual use, the temperature may vary due to different ambient temperatures, which is a normal phenomenon.
Product Name	Recovery Therm Cube	Product Model	HXR-A130																			
Product size	149mm*76mm*50mm	Rated Power	20W																			
Battery capacity	3200mAh	Charging Voltage	5V																			
Product Weight	387g	Charging Port	USB Type-C																			

Applicable body parts		Wearing method	Product charging	Troubleshooting	Maintenance and storage methods													
 Head  Neck  Waist  Thighs  Shoulders  Arms  Knees  Calves		Step 1: Secure the strap to one side of the product with the velcro. (Figure 1) Step 2: Bend the product and fit it onto the application area. (Figure 2) Step 3: Wrap the strap around the application area and fix it to the other side of the product with the velcro. (Figure 3)    ①②③ <i>*Note:</i> When used on the neck, shoulders or knees: It can be bent as needed and then attached to the application area.	When used on the arms or legs, the “strap frame” module can be installed first, and the user experience can be enhanced by changing the fixed position of the strap. The installation process is as shown in the following figure: Step 1: Install the strap frame on the product and fix the strap to one side of the product strap frame the velcro (Figure 1) Step 2: Attach the product to the application area and wrap the strap around it. (Figure 2) Step 3: Fix the velcro to the other side of the product. (Figure 3)    ①②③	If it is used for the first time or has not been used for a long time, please charge it as follows: - The USB cable in the accessory should be plugged into the charging socket or charging port which is available on the bottom left side of the device - Insert the USB cable into the power adaptor and connect the power adaptor to the electricity  - After the charging starts, the power indicator flashes on the device - The device will be fully charged for about 3.5 hours, continue charging for 3 seconds after the power indicator goes off, please unplug the power in time after the full charging of the device. <i>*Note:</i> Do not use this device while charging.	<table><tr><th>Possible conditions</th><th>Possible causes/handling methods</th></tr><tr><td>Cannot start after pressing the power button</td><td>Please check whether the product has sufficient power or not Please confirm whether the product start-up operation is standardized (press and hold any function button for more than 1 second)</td></tr><tr><td>Stopped during its use</td><td>The device will automatically stop working after 15min of use, which is a normal phenomenon, if you need to continue to use the power can be turned on again Please confirm whether the device is in low battery status, if the battery power is insufficient, please charge it and use it again after its fully charged.</td></tr><tr><td>Unable to charge or very short life span</td><td>If the ambient temperature is particularly high or low, the device will not start, please charge at 10-40 °C. Please check whether the USB connection cable and the host power adapter are in good contact The first charge or long time no use, please first full charge before its use. The device cannot be charged during startup, please turn off the power before charging</td></tr><tr><td>Redness of the skin after use</td><td>This is a normal phenomenon after the skin has been cold/heated, and will return to normal within half an hour after you stop using the device. Please check if your skin is allergic to the material of the massage head</td></tr><tr><td>Device heating up after shutdown</td><td>This is the normal phenomenon after the end of the usage of the device in the cool function, and can be restored to normal temperature after a few minutes of stopping the device.</td></tr></table> <i>*Note:</i> If the problem still exists after following the above steps, please cut off the power and contact customer service to deal with it.	Possible conditions	Possible causes/handling methods	Cannot start after pressing the power button	Please check whether the product has sufficient power or not Please confirm whether the product start-up operation is standardized (press and hold any function button for more than 1 second)	Stopped during its use	The device will automatically stop working after 15min of use, which is a normal phenomenon, if you need to continue to use the power can be turned on again Please confirm whether the device is in low battery status, if the battery power is insufficient, please charge it and use it again after its fully charged.	Unable to charge or very short life span	If the ambient temperature is particularly high or low, the device will not start, please charge at 10-40 °C. Please check whether the USB connection cable and the host power adapter are in good contact The first charge or long time no use, please first full charge before its use. The device cannot be charged during startup, please turn off the power before charging	Redness of the skin after use	This is a normal phenomenon after the skin has been cold/heated, and will return to normal within half an hour after you stop using the device. Please check if your skin is allergic to the material of the massage head	Device heating up after shutdown	This is the normal phenomenon after the end of the usage of the device in the cool function, and can be restored to normal temperature after a few minutes of stopping the device.	● Maintenance method - Be sure to unplug the power cable or cord before maintenance. - Do not use detergent, thinner, gasoline or other chemicals to wipe. - When the device is dirty, it can be wiped with a dry soft cloth or dry towel to avoid washing in running water or other liquids. - After each use, please clean the massage head in time. When the massage head is particularly dirty, wipe or disinfect it with a soft cloth dipped in a small quantity of medically - approved cleaner. - Please do not disassemble, repair or modify this product without permission to avoid accidents or failures. ● Storage method - Be sure to unplug the power supply during storage. - Please do not keep it in a damp place or a place exposed to direct sunlight. - If the product is not used for a long time, please fully charge the product, put it in the storage bag, and then place it in a dry and ventilated place. - Do not store it in a place where infants or pets can touch it. - Please do not place it on cloth, paper products or wooden furniture that are easy to fade.
Possible conditions	Possible causes/handling methods																	
Cannot start after pressing the power button	Please check whether the product has sufficient power or not Please confirm whether the product start-up operation is standardized (press and hold any function button for more than 1 second)																	
Stopped during its use	The device will automatically stop working after 15min of use, which is a normal phenomenon, if you need to continue to use the power can be turned on again Please confirm whether the device is in low battery status, if the battery power is insufficient, please charge it and use it again after its fully charged.																	
Unable to charge or very short life span	If the ambient temperature is particularly high or low, the device will not start, please charge at 10-40 °C. Please check whether the USB connection cable and the host power adapter are in good contact The first charge or long time no use, please first full charge before its use. The device cannot be charged during startup, please turn off the power before charging																	
Redness of the skin after use	This is a normal phenomenon after the skin has been cold/heated, and will return to normal within half an hour after you stop using the device. Please check if your skin is allergic to the material of the massage head																	
Device heating up after shutdown	This is the normal phenomenon after the end of the usage of the device in the cool function, and can be restored to normal temperature after a few minutes of stopping the device.																	