

1. Reclaim

At the end of the service life of the product, do not discard it together with general household waste. The product should be placed in the recycling bin designated by the government, or the product can be sent back to the after-sales service center of the original factory. The personnel of the service center will deal with it in an environmentally sound manner.

1. Maintenance, cleaning, storage and transportation

01. Clean the equipment with a slightly wet towel and dry it with a dry soft cloth.

02. If the equipment is not used for a long time, please store it in the packing box.

03. When the equipment is not used for a long time, the battery shall be charged once a month.


04. The massage instrument shall be stored in a well ventilated warehouse without corrosive gas, 250mm above the ground;

05. It is strictly prohibited to stack massage instruments in more than 5 layers.

06. When the massage instrument is not used for a long time, it should be packaged and stored after cleaning. It is best to take it out and charge it once a month to prevent moisture.

07. This product is applicable to automobile, train, ship and aircraft transportation;

08. During dynamic transmission, the products shall be fixed to prevent violent vibration, collision, tumbling and falling, and shall be upward, handled with care and rainproof;

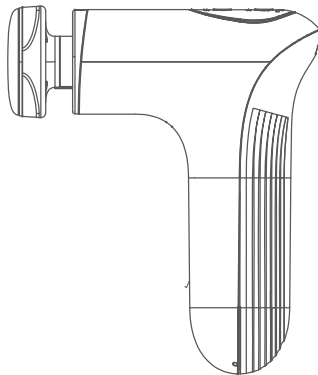

Mini Thermal Massage Gun

HXR-B029 User's manual


Thank you for purchasing our products. Before use, please read this manual carefully to ensure correct usage. Please keep this manual properly for future query.

1- Introduction and application of massage gun:

- Thermal massage gun is a device for deep relaxation of muscle fascia, which is suitable for physical massage after exercise and fitness. It can reduce muscle stiffness and pain, relieve muscle fatigue and accelerate blood circulation, so as to improve the health of soft tissue.

Purpose

- Suitable for physical massage and skin spa beauty care before / after exercise and fitness

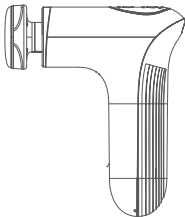

Thermal massage gun ---- your health care master

Important information

In order to reduce the risk of failure, electric shock, injury, fire, death and equipment damage, please carefully read all instructions and warnings about the product in this manual before using the product.



Prohibition

1. It is forbidden for minors to use. If they are injured, it is forbidden to use this sample. Consult your doctor before using this device in case of muscle injury.
2. It is forbidden to use on the surface of non dry and unclean skin. It can be used across clothes. Gently press and move each part on the skin for about 60 seconds.
3. Do not immerse in water or let water enter the ventilation hole of the product.
4. It is forbidden to place any object in the ventilation hole of the product.
5. It is forbidden to put the product into the fire or close to the fire source.
6. It is forbidden to drop the product intentionally.
7. It is forbidden to change or modify the product in any way.
8. It is forbidden to operate or charge the product unattended.
9. Fingers, hair or other parts of the body shall not be close to the shaft and back of the massage head to avoid pinch injury.
10. It is prohibited for children under 8 years old. Children shall not use, clean and maintain the product without supervision.
11. The equipment has a heating surface, so people who are not sensitive to heat should use it carefully.



illustration: Do not touch this part of the mold when starting the machine after installing the massage head to avoid clamping injury.



illustration: When installing the massage head and starting the machine, do not touch this part of the film to avoid pinch injury.

2. Declaration

Consult your doctor before use if you have any of the following conditions.


1. Pregnant mothers and cancer patients
2. Heart disease, cardiac pacemaker or defibrillator
3. Recent surgical treatment
4. Aneurysm, bleeding, use of blood thinner
5. There are differences in sensitive areas, head or joint bones
6. Complications caused by other diseases (such as neuropathy, etc.)

3. Precautions for use

01. After charging, unplug the charger.
02. It is forbidden to use this product under charging state.
03. Please use the product correctly according to the instructions.

1.Product introduction


Thank you for choosing e-cozy Thermal Massage Gun

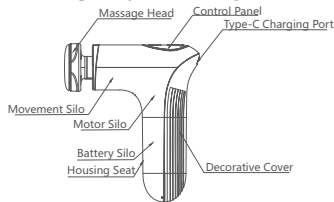

Massage Head Introduction

- The **Heating massage head** is mainly used for relaxing the whole body and shaping skin beauty.
- The **Spherical head** is mainly used for large muscle groups.
- Cylindrical head** is mainly used to impact deep tissue.
- The **U-shaped head** is mainly used on both sides of the neck and spine.
- Flat head** is mainly used for latissimus dorsi and trapezius. Usage operation

2- Application method

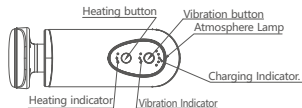

1. Long press the vibration mode key for 2 seconds, and the vibration function is turned on. At this time, the force is low, and the indicator light turns blue; Short press once, the force is medium, and the indicator light turns purple; Short press again, the strength is high-grade, and the indicator light turns red; Short press again to shut down or long press the key for 2 seconds to shut down.

2. Long press the heat button for 2 seconds, the heating function is turned on, the temperature is about 35 °C, and the indicator light turns blue; Short press once, the temperature is about 40 °C, and the indicator light turns purple; Short press again, the temperature is about 45 °C, and the indicator light turns red; Short press the close thermal key again or long press the thermal key for 2 seconds to close.

3. To have the heating function of point 2, heating massage head must be installed, and other massage heads do not have this function.

4. Vibration and heating functions can be used separately.

5. The automatic timing of vibration and heating functions is 15 minutes.

Assembling Massage Head




Note: The black solid line is the product alignment line, the actual product is bright black here

Installation of massage head:

1. The disassembly and assembly of massage head are in-line, without screwing.
2. The red straight bar of cold and hot head is aligned with the gap of piston.
3. Be sure to plug it to the bottom

Figure for installation of heat head

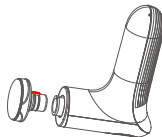


Error Reporting Notes:

Charging:

The charging port is located at the top of the product

Handle:

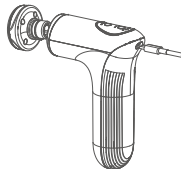

1. Before initial use, charge for about 4 hours and fully charge.
2. For the charging battery, a special USB cable will be provided to connect the device type-C socket, and any 5-volt charger head will be used at the other end. (5V, 2A recommended)
3. A II power indicator lights on the instrument flicker regularly, indicating that they are charging.
4. When all power indicator lights on the instrument stop flashing and remain on, it indicates that the power is fully charged.
5. The battery can be charged at any time and at any power level under the specified environmental conditions.
6. When the power is lower than 25%, it shall be charged in time.
7. It is not recommended to completely drain the battery to 25%.
8. The average running time of the equipment after loading is three hours, which depends on the vibration frequency and applied pressure during use.

Basic Parameters

Product Name: Mini thermal massage gun	Product model: HXR-B029
Net weight: 457g	Input voltage: 5V, 2A
Power: 10W (maximum 15W)	Battery Specification: 7.4V 2200mAh
Speed gear: 3rd gear adjustment	Charging time: about 4 hours
Service environment: 0-40 °C	Product size: 151 × 48 × 158MM
Battery life: about 4-5 hours Heating for about 5-6 hours	Automatic timing: knock for 15 minutes Heating for 15 minutes

Temperature value:

Heating: 35 degrees for 1st gear, 40 degrees for 2nd gear and 45 degrees for 3rd gear

The above deviation will be plus or minus 5 degrees depending on the environment and use conditions.

Product executive standard: GB4706.1-2005 GB4706.10-2008

3- Common problem

Unable to power on:

1. Over heating

When the fascia gun is used in an environment above 40 °C, the rib mold gun may trigger the temperature overheating protection, and the host cannot start normally.

Solution: stand still for 2 hours in an environment below 40 °C, and then turn it on for use.

2. Battery level is "0":

When the battery is out of power, if the user does not charge for a long time, the battery will produce low-power protection, collectively referred to as over discharge protection. At this time, the fascia gun will not work normally.

Solution: keep the power off mode, connect the charger, let the fascia gun enter the charging mode, unplug the charger after being fully charged, and then start the fascia gun to work normally.

Unable to charge:

1. Massage gun temperature too high:

When the temperature of the main body of the massage gun is too high, it will not be charged when connected to the charger.

Solution: Stand still for 2 hours in an environment below 40 °C, and then charge.

2. Charger temperature too high

When massage gun is charged in an environment above 40 °C, the charger may trigger overheating protection, and the massage gun cannot be charged normally.

Solution: After standing still in the environment below 40 °C for 2 hours, then charge it.

Troubleshooting of Gun operation

Fault description	Troubleshooting
1. The massage head slides out and falls	Check whether the massage head is installed in place (see page 7 of the manual for specific operation)
2. The product cannot be started	Check whether the power is sufficient (if it has not been used for a long time, please charge it for more than 6 hours)
3. The product is unable to get charge	Check whether the charging port is inserted in place
4. Suddenly stop after continuous operation for 10 minutes	This phenomenon is the setting of long-time operation protection function. It is normal to change from operation to standby mode and restart by pressing the key for a short time
5. The knocking function stops after the massage head is blocked	This phenomenon is the setting of abnormal operation protection function. It is normal to change from operation to standby mode and restart by pressing the key for a short time