

B260 英文通用说明书, 折页尺寸: 70*110mm

HANDHELD HEATING MASSAGER

Instruction Manual



Safety instructions

Before using this product, please be sure to read and follow the operating instructions in the "Operating Instructions" carefully. Violation of the instructions may result in danger or damage to the equipment.

Be sure to unplug the product from the outlet before using or cleaning it. Please do not throw the product into water or other liquids.

If the product accidentally falls into the water while charging, please pull out the plug immediately and do not touch the product directly into the water.

This product contains lithium batteries, please do not disassemble, modify or throw in the fire to avoid danger.

Caution

This product has a heating function, the heat source is not sensitive or easy to heat allergic users, please use the heating function with caution.

This product has a heating function, please do not repeatedly and continuously use in the same part to avoid burns.

Please do not plug and unplug the USB charging cable with wet hands to avoid electric shock or accidents.

This product is not a substitute for drugs and medical devices.

Please do not place the product in front of an unattended child.

Children or those with limited ability should use this product under adultsupervision.

Please stop using the product immediately if you feel uncomfortable during use.


Product parameters

Product Name	Handheld Heating Massager	Product Model	HXR-B260
Product size	200mm*130mm*47mm	Rated Power	8 W
Battery capacity	3.7V 1800mAh	Charging Voltage	5 V
Product Weight	196g	Charging Port	USB Type-C


Product list



Main Device* 1

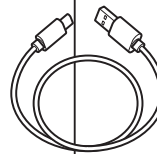


Type-C Charging Cable* 1

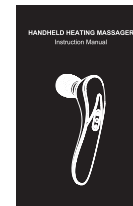

Instructions* 1

Device description

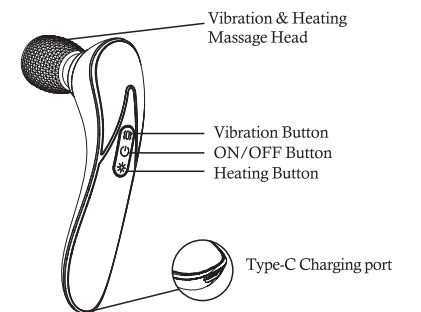

Usage instructions

Long press the on/off button for 2 seconds to turn on the massager, the button indicator lights up and it is under the first level of vibration /heating. The massager's single working time is 15 minutes. After 15 minutes, it automatically turns off. Long press the on/off button or vibration button for 2 seconds under the power-on state can also turn off the massager.

Under power-on state, short press vibration button to adjust the intensity, total 3 levels.

Level 1, button indicator light is blue.


Level 2, button indicator light is purple.

Level 3, button indicator light is red.

Device charging

If it is going to be used for the first time or has not been used for a long time, please charge it as per given instructions:

1. The USB cable in the accessory should be plugged into the charging socket on the bottom left side of the product
2. Insert the USB cable into the power adapter and connect the power.

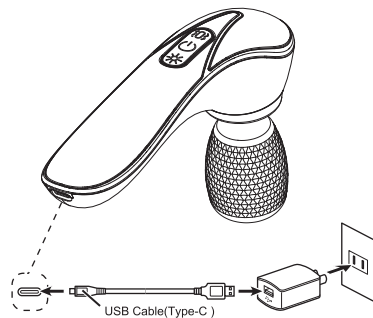

3. During charging, the indicator light of the ON/OFF keeps red blinking.
4. It takes about 2-3 hours for the massager to be fully charged, and the indicator light of the on/off button turns green always on after it is fully charged, please unplug the power supply in time after it is fully charged.

*Note: Do not use this device during charging.

Troubleshooting

Possible conditions	Possible causes/handling methods
Unable to start after pressing the power key	Please confirm whether the power of the device is sufficient.
	Please confirm whether the massager startup operation is normal (long press the on/off button for 2 seconds)
Equipment not heating/vibration	Please make sure that the hot massage head is properly pressed on the skin surface
	Please confirm whether the power of the device is sufficient.
Stopped during its use	The device will automatically stop working 15 minutes after use, which is a normal phenomenon. If you need to continue using it, you can turn on the power again.
	Please confirm whether the device is in a low power state. If the battery is low, please charge it before use.
Unable to charge or short battery life	If the device cannot be started at a particularly high or low ambient temperature, please charge it at 10-40°C.
	Please check whether the USB connection cable and the host power adapter are in good contact.
	When charging for the first time or not using for a long time, please fully charge before using.

*Note: If the problem still exists after following the above steps, please cut off the power and contact the customer service.

Maintenance and storage method

1. Before maintenance, please make sure to unplug the power cord.
2. Do not use detergent, thinner, gasoline and other chemicals to wipe.
3. When the host is dirty, use a dry soft cloth or dry towel to wipe, avoid cleaning in running water or other liquids.
4. Please do not disassemble, repair or modify this product by yourself to avoid accidents or malfunction.
5. Please do not store in a humid place or a place with direct sunlight.
6. When the product is not used for a long time, please put the product in a dry and ventilated place after it is fully charged.
7. Please do not store in places where babies and children, pets can touch.