



折页，双面，单页尺寸：105*50MM

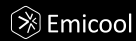

INSECT BITE HEALER

Instruction Manual

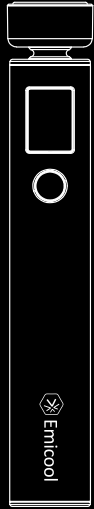

 Safety instructions

Before using this product, please be sure to read and follow the operating instructions in the "Operating Instructions" carefully. Violation of the instructions may result in danger or damage to the equipment.

Be sure to unplug the product from the outlet before using or cleaning it. Please do not throw the product into water or other liquids.

If the product accidentally falls into the water while charging, please pull out the plug immediately and do not touch the product directly into the water.

This product contains lithium batteries, please do not disassemble, modify or throw in the fire to avoid danger.



 Caution

This product has a heating function, the heat source is not sensitive or easy to heat allergic users, please use the heating function with caution.

This product has a heating function, please do not repeatedly and continuously use in the same part to avoid burns.

Please do not plug and unplug the USB charging cable with wet hands to avoid electric shock or accidents.

This product is not a substitute for drugs and medical devices.

Please do not place the product in front of an unattended child.


Children or those with limited ability should use this product under adult supervision.

Please stop using the product immediately if you feel uncomfortable during use.


Product Parameters

Product Name	INSECT BITE HEALER	Product Model	HXR-B270
Product size	20mm*20mm*118mm	Rated Power	3 W
Battery capacity	3.7V 260mAh	Charging Voltage	5 V
Product Weight	30g	Charging Port	USB Type-C

Product List

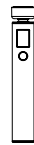

Main Device*1

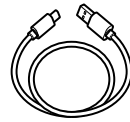

Type-C Charging Cable*1




Lanyard* 1



Instructions* 1

Device description

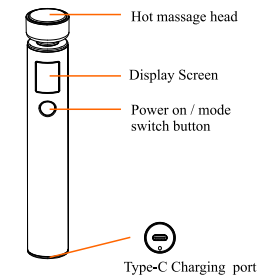

Proper use

Intended use

The insect bite healer is intended for relieving itching and swelling caused by the insect bites or stings by applying heat (local hyperthermia).

Target group

Only use the insect bite healer to treat insect bites or stings. The insect bite healer is intended for domestic/private use only, not for commercial use.

Indication/clinical benefits

Many insect bites and stings result in painful and swollen sting and bite marks. These bites/stings often cause severe itching. Itching and swelling can be alleviated by the targeted effect of heat on the affected area.

Usage instructions

Short press the power on button, the display shows A and enter the 1st gear, continue to short press the power on button to switch gears, the itching stopper has a total of 3 gears:

- 1st Mode: the screen display A, temperature 45°C, vibration on, working time one minute;
- 2nd Mode: the screen display B, temperature 50°C, vibration on, working time one minute;
- 3rd Mode: the screen display C, temperature 55°C, vibration on, working time one minute;

*Note: The default single working time of the anti-itching instrument is one minute, and it will automatically shut down after one minute.

Use the following method:

1. Open the insect bite healer, select the desired mode
2. Press the hot massage head on the mosquito bite to hold for one minute
3. When the insect bite healer stops working, take it off

Warnings and safety notes


- Do not use the device if you suffer from diabetes. Sensitivity to pain could, in certain circumstances, be reduced and diabetes may burn themselves.
- Do not use the device on people with disabilities, children under 3 years old or people with reduced sensitivity to heat. This includes people with skin alterations due to illness or scarred tissue in the application area, after taking pain relief medication, drinking alcohol or taking drugs. Do not use the device on acutely or chronically diseased (injured or irritated) skin (e.g. inflamed skin, whether painful or not, reddened skin, rashes, (e.g. allergies), burns, frostbite, bruises, swellings, both open and healing wounds, and post-operative scars where the healing process could be affected).
- Do not use the device on sensitive areas of skin.
- Do not use the device if you have any sensory impairment that reduces the feeling of pain (e.g. metabolic disorders).
- Do not use the device if you are using lotions, creams and gels locally at the same time.
- Do not use the device if you suffer from persistently irritated skin due to long-term heat application on the same area of skin.

Device charging

If it is going to be used for the first time or has not been used for a long time, please charge it as per given instructions:

1. The USB cable in the accessory should be plugged into

1. The USB cable in the accessory should be plugged into the charging socket on the bottom left side of the product
2. Insert the USB cable into the power adapter and connect the power.

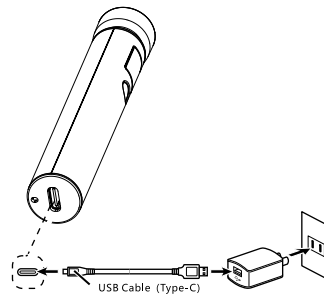

3. When charging starts, the white light above the product button flashes
4. It takes about 90 minutes to fully charge the device, and the four lights above the button of the itching stopper are always on after full charging, and will go off after 3 minutes of continued charging, please unplug the power in time after full charging the device.

*Note: Do not use this device during charging.

Troubleshooting

Possible conditions	Possible causes/handling methods
Unable to start after pressing the power key	Please confirm whether the power of the device is sufficient.
	Please make sure that the hot massage head is properly pressed on the skin surface.
Equipment not heating/vibration	Please confirm whether the power of the anti-itch instrument is sufficient.
	The device will automatically stop working one minute after use, which is normal phenomenon. If you need to continue using it, you can turn on the power again.
Stopped during its use	Please confirm whether the device is in a low power state. If the battery is low, please charge it before use.
	If the device cannot be started at a particularly high or low ambient temperature, please charge it at 10-40 °C.
Unable to charge or short battery life	Please check whether the USB connection cable and the host power adapter are in good contact.
	Please check whether the USB connection cable and the host power adapter are in good contact.
	When charging for the first time or not using for a long time, please fully charge before using.

*Note: If the problem still exists after following the above steps, please cut off the power and contact the customer service.

Maintenance and storage method

1. Before maintenance, please make sure to unplug the power cord.
2. Do not use detergent, thinner, gasoline and other chemicals to wipe.
3. When the host is dirty, use a dry soft cloth or dry towel to wipe, avoid cleaning in running water or other liquids.
4. Please do not disassemble, repair or modify this product by yourself to avoid accidents or malfunction.
5. Please do not store in a humid place or a place with direct sunlight.
6. When the product is not used for a long time, please put the product in a dry and ventilated place after it is fully charged.
7. Please do not store in places where babies and children or pets can touch.