


DEEP TISSUE VIBRATION MASSAGER

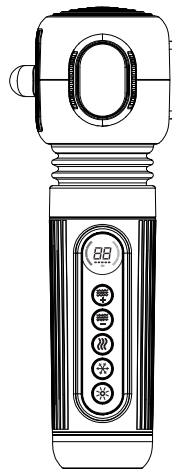

B510 USER'S MANUAL

Safety instructions

Before using this product, please be sure to read and follow the operating instructions in the "Operating Instructions" carefully. Violation of the instructions may result in danger or damage to the equipment.

Be sure to unplug the product from the outlet before using or cleaning it.

Please do not throw the product into water or other liquids.

If the product accidentally falls into the water while charging, please pull out the plug immediately and do not touch the product directly into the water.

This product contains lithium batteries, please do not disassemble, modify or throw in the fire to avoid danger.

Caution

This product has a heating function, the heat source is not sensitive or easy to heat allergic users, please use the heating function with caution.

This product has a heating function, please do not repeatedly and continuously use in the same part to avoid burns.

Please do not plug and unplug the USB charging cable with wet hands to avoid electric shock or accidents.

This product is not a substitute for drugs and medical devices.

Please do not place the product in front of an unattended child.


Children or those with limited ability should use this product under adultsupervision.

Please stop using the product immediately if you feel uncomfortable during use.

Product parameters

Product Name	Deep Tissue Vibration Massager	Product Model	HXR-B510
Product Size	216mm*78mm*91mm	Rated Power	20W
Battery Capacity	7.4V 2200mAh	Charging Voltage	5V $\overline{\text{---}}$ 2A
Product Weight	600g	Charging Port	USB Type-C

Product list

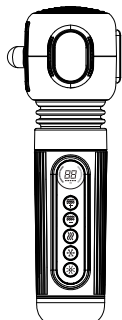

Main Device* 1




Type-C Charging Cable* 1





Instructions* 1

Device description

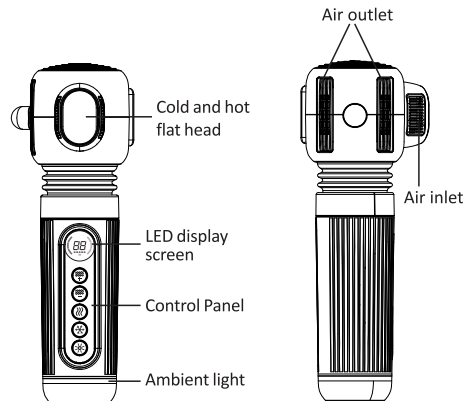

Silicone massage head

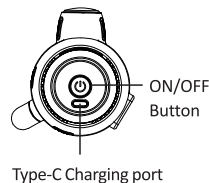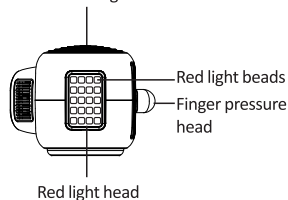

Massage head function

Finger pressure head


Point pressing is suitable for deep tissue or concentrated pain point massage, such as the palms, soles of the feet, and deep muscles of the calves.



(X signifies the area for release, red indicates the area of pain)

Cold and hot flat head

It is suitable for all muscle groups throughout the body, such as large areas like the waist and back, buttocks, and the front of the thighs.

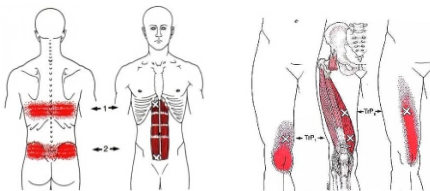

(X signifies the area for release, red indicates the area of pain)

Red light head


It is applicable to scenarios such as joint repair, sports recovery, and chronic inflammation management.



(X signifies the area for release, red indicates the area of pain)

Silicone massage head


Made of silicone material, it is suitable for sensitive areas or areas with thinner muscles (such as pectoral muscles, shoulders and neck), with gentle force and less likely to cause irritation.

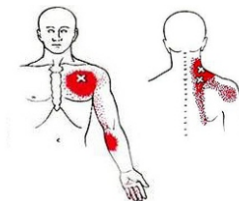

(X signifies the area for release, red indicates the area of pain)

Device operation instructions

Control Panel

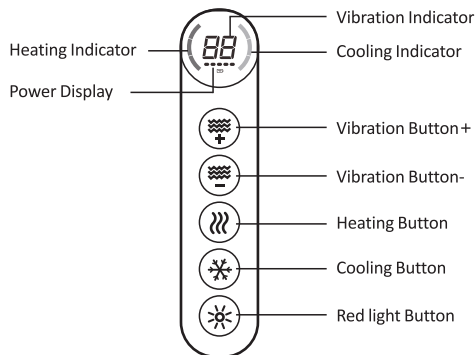

Press the power button briefly, and the ambient light at the bottom will light up. At this time, the massage device is in standby mode. If there is no operation within one minute, it will automatically shut down.

When the massage device is in standby mode, the functions can be adjusted through the function keys. The functions of each key are as follows:

 **Vibration Button+** : Short press can increase the vibration intensity, up to 10 levels at most

 **Vibration Button-**: Short press can decrease the vibration intensity, and the lowest level can be reduced to 0

 **Heating Button** : Short press can switch between the hot compress Settings. There are a total of three hot compress Settings, with temperatures ranging from 40°C-45°C-50°C

 **Cooling Button** : Short press can switch between the cold compress Settings. There are a total of three cold compress Settings, with temperatures ranging from 20°C-15°C-10°C

 **Red light Button** : Short press can turn on/off the red light mode. The red light consists of 660nm light and 850nm infrared light


The massage device will automatically turn off after being turned on for 15 minutes. When it is powered on, press and hold the power button for 2 seconds to turn off the massage device. If you need to continue using it, please turn it on again.

*Note: The hot/cold compress function is affected by the ambient temperature, and the actual temperature may have a deviation of $\pm 3^{\circ}\text{C}$.
When the device is in cold compress mode, please ensure that the air outlet and air inlet are not blocked to avoid , affect the cooling effect of the equipment.

Device charging

If it is going to be used for the first time or has not been used for a long time, please charge it as per given instructions:

1. The USB cable in the accessory should be plugged into the charging socket on the bottom left side of the product
2. Insert the USB cable into the power adapter and connect the power.

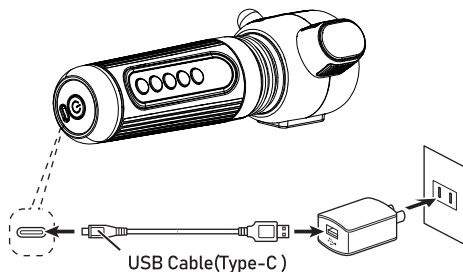

3. After the charging starts, the battery level indicator light on the LED screen of the massage device flashes.
4. It takes approximately 3 hours to fully charge the product, after full charging, the LED battery level display of the massage device will show 5 bars, please unplug the power supply in time after full charging.

*Note: Do not use this device during charging.

Troubleshooting

Possible conditions	Possible causes/handling methods
Unable to start after pressing the power key	Please confirm whether the power of the device is sufficient.
	Please confirm whether the massager startup operation is normal(Short press the on/off button)
Equipment not heating/vibration	Please make sure that the hot massage head is properly pressed on the skin surface
	Please confirm whether the power of the device is sufficient.
Stopped during its use	The device will automatically stop working 15 minutes after use, which is a normal phenomenon. If you need to continue using it, you can turn on the power again.
	Please confirm whether the device is in a low power state. If the battery is low, please charge it before use.
Unable to charge or short battery life	If the device cannot be started at a particularly high or low ambient temperature, please charge it at 10-40℃.
	Please check whether the USB connection cable and the host power adapter are in good contact.
	When charging for the first time or not using for a long time, please fully charge before using.

*Note: If the problem still exists after following the above steps, please cut off the power and contact the customer service.

Maintenance and storage method

1. Before maintenance, please make sure to unplug the power cord.
2. Do not use detergent, thinner, gasoline and other chemicals to wipe.
3. When the host is dirty, use a dry soft cloth or dry towel to wipe, avoid cleaning in running water or other liquids.
4. Please do not disassemble, repair or modify this product by yourself to avoid accidents or malfunction.
5. Please do not store in a humid place or a place with direct sunlight.
6. When the product is not used for a long time, please put the product in a dry and ventilated place after it is fully charged.
7. Please do not store in places where babies and children, pets can touch.