

B610英文说明书

铜版纸彩色印刷

单页尺寸：120mm*85mm



BELLY MASSAGER

INSTRUCTIONAL MANUAL



HXR-B610

Danger

Before using this HXR - B610 Belly Massager, please carefully read and strictly comply with the operating instructions in this manual. The manufacturer shall not be liable for any consequences arising from incorrect use or failure to follow these instructions.

Do not use this product on open wounds, as it may cause infection and aggravate the injury.

Avoid using this massager in combination with other electronic medical devices simultaneously without professional guidance, as it may cause electromagnetic interference and affect the normal operation of medical devices.

This equipment should be used within the temperature range of 0°C - 40°C. Store the equipment and its accessories within the temperature range of - 20°C - 60°C. If the ambient temperature is too high or too low, it may cause equipment malfunction or damage.

Caution

- Do not use this product while driving or operating machinery to avoid distractions that may lead to accidents.
- Keep this massager away from water sources (such as bathtubs, sinks) and open flames to prevent electrical short - circuits or fire hazards.
- Do not use this product in environments with strong electromagnetic interference (such as near high - power electrical equipment), high temperature (such as direct sunlight in summer), or high humidity (such as bathrooms without proper ventilation), as it may affect the performance and service life of the product.
- Do not attempt to disassemble, repair, or modify this product by yourself. If the device fails to start or malfunctions, please contact the after - sales service in time.
- When plugging or unplugging the Type - C charging cable, make sure your hands are dry to avoid the risk of electric shock.
- This product is not suitable for patients with severe skin diseases (such as severe eczema, psoriasis), acute abdominal pain, or abdominal tumors.
- Do not use it for patients who have recently undergone abdominal surgery, as it may affect wound healing.
- Do not use it for patients with symptoms such as abdominal aortic aneurysm, as improper use may lead to serious consequences.
- Before use, ensure that there are no metal objects, such as belts with large metal buckles, in the area to be massaged, and remove accessories such as body jewelry that may interfere.
- This is a personal care massage product and cannot replace drugs and medical devices. If you have skin allergies, skin diseases, or other health problems, please use this product under the guidance of a doctor.
- If you feel any discomfort during use, please stop using it immediately.
- Please do not use this product for purposes other than personal care massage as described in the manual to avoid accidents and product damage.
- This product contains a lithium - ion battery. If the battery needs to be replaced, please contact the official after - sales service. Do not disassemble or replace the battery by yourself.
- This product is made of plastic materials. Please do not discard it at will. Do not throw it into fire to avoid danger.
- When discarding this product, please handle it in accordance with local environmental protection regulations to avoid environmental pollution.

Product parameters

Device Name	Belly Massager	Device Model	HXR-B610
Device Size	151mm*141mm*138mm	Rated Power	10W
Battery Capacity	7.4V 1500mAh	Charging Voltage	5V, 2A
Device Weight	663g	Charging Port	USB Type-C


*Note: Because the belt has good elasticity, so there may be ±5CM error between the measurement size and the actual size.

What is included in the package ?




Main Device*1 Type-C Charging Cable*1 Instructions Manual*1

Interchangeable Massage Heads*4

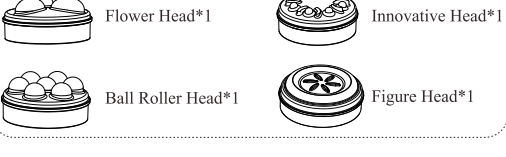

Functions and names of each part

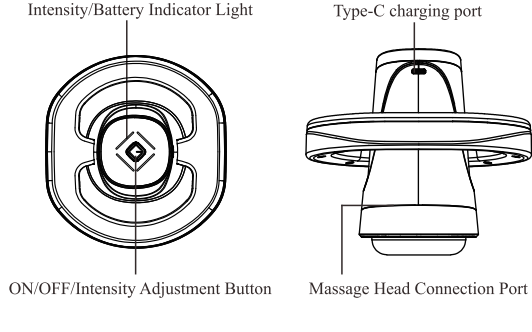

Main Unit Parts

Intensity/Battery Indicator Light:
The white indicator light is used to display the remaining battery level or the current intensity.

ON/OFF/Intensity Adjustment Button:
Powers the device on or off.
Short - press to adjust the 4 different vibration intensities.

Type - C Charging Port:
For connecting the charging cable to recharge the 1500mAh battery.

Massage Head Connection Port:
Where the interchangeable massage heads are attached.

Massage Head Installation

- Flower Head:**
Designed for a gentle, all - around massage experience, suitable for areas like the abdomen and thighs to help with skin tightening and reducing stretch marks.
- Ball Roller Head:**
Features multiple balls for a rolling massage effect, ideal for relieving muscle soreness in larger muscle groups such as the shoulders and calves.
- Innovative Head:**
Offers a unique massage pattern, great for body sculpting and cellulite massage on areas like the waist and arms.
- Figure Head:**
Shaped to fit specific body contours, perfect for targeted massage on smaller areas or for more precise pressure application.

Product operating instructions

Operation Instructions

Power On: Press and hold the power key for 2 seconds to turn on the device. It will default to Speed 1, with 1 white light on.

Switch Speeds: When the device is on, short - press the power key to cycle through the 4 speeds. Speed 1 has 1 white light on, Speed 2 has 2 white lights on, and so on (Speed 4 has 4 white lights on).

Auto - Shutoff: The device will automatically turn off after 15 minutes of continuous operation.

Power Off: When the device is working, press and hold the power key for 2 seconds to turn it off.


Check Battery (When Off): With the device off, short - press the power key to check the battery:

- 1 light flashing: Battery <5%
- 1 light on: Battery 6% - 25%
- 2 lights on: Battery 26% - 50%
- 3 lights on: Battery 51% - 75%
- 4 lights on: Battery 76% - 100%

Charging Instructions

During charging, the battery status is shown by white lights:
Battery 0 - 25%: 1st white light flashes
Battery 26 - 50%: 1st white light on, 2nd white light flashes
Battery 51 - 75%: 1st & 2nd white lights on, 3rd white light flashes
Battery 76 - 99%: 1st, 2nd & 3rd white lights on, 4th white light flashes
Battery 100%: All 4 white lights on, then turn off automatically after 5 minutes.

Attaching/Replacing message heads

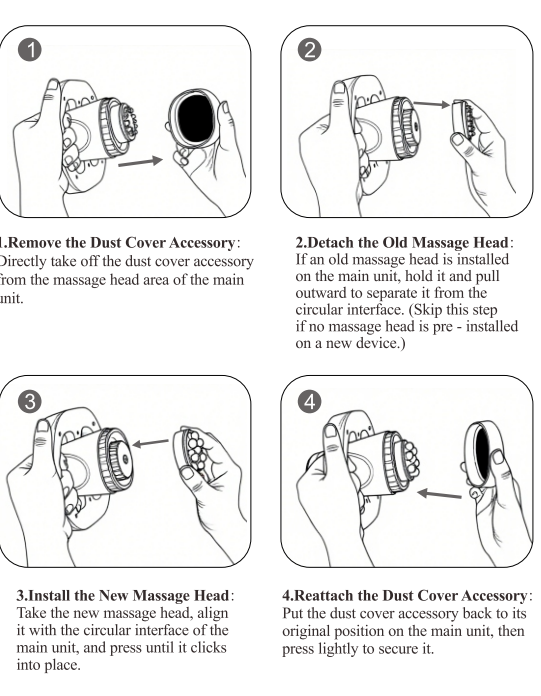

Wearing method



Troubleshooting

Possible Conditions	Possible Causes/Handling Methods
Device Not Powering On	1. Check the Battery: Ensure the device is charged. Connect the Type - C charging cable to a 5V power source and let it charge for at least 30 minutes. If it was completely dead, it may take a few minutes for the power indicator (if any) to show signs of life. 2. Check the ON/OFF Button: Make sure you are pressing the button correctly (press and hold for power - on/off as per instructions). Sometimes, dirt or debris around the button may interfere; gently clean the area around the button.
Massage Head Not Attaching Properly	1. Alignment Issue: Ensure the massage head is properly aligned with the connection port. Check for any debris or damage on the port or the massage head that may be preventing a secure attachment. 2. Mesh Cover Interference: Make sure the mesh cover is removed before attempting to attach or remove the massage head. A misaligned or improperly placed mesh cover can also affect the attachment.
Vibration Intensity Not Changing	1. Button Function: Check if you are short - pressing the Vibration Level Button correctly. If the button seems unresponsive, try powering the device off and on again to reset it. 2. Battery Level: A very low battery may affect the functionality of the intensity adjustment. Charge the device fully and then test the intensity adjustment again.

*Note: If the problem still exists after following the above steps, please cut off the power and contact customer service to deal with it.

Maintenance and storage methods

● Maintenance method

Cleaning the Main Unit: Use a soft, dry cloth to wipe the surface of the main unit. Do not use water or any liquid cleaners directly on the main unit to avoid damage to the electrical components.

Cleaning the Massage Heads: Remove the massage heads and clean them with a mild soap and water solution (if they are not electrical components). Rinse them thoroughly and dry them completely before reattaching to the main unit. Make sure to also clean the mesh covers in the same way.

● Storage method

Charging: Before long - term storage, charge the device to about 50 - 70% battery level. Do not store it with a fully depleted battery.

Environment: Store the massager in a cool, dry place. Avoid direct sunlight, high temperatures, and high humidity environments, as these can damage the battery and other components.

Disassembly for Storage: If desired, you can remove the massage heads and store them separately in a small, clean container to prevent damage and keep them organized.