


三折页双面印刷

单页尺寸：80\*160MM

FOOT MASSAGER

HXR-C060 USER'S MANUAL

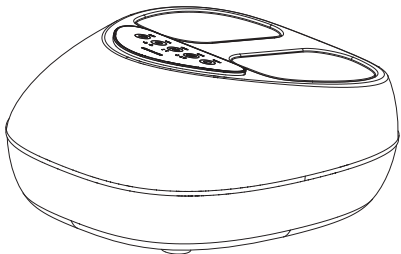

Safety guidelines

1.Description of important matters

- In order to ensure safe and correct use of the product, please be sure to check carefully before use Operation Manual.
- Be sure to use this product in accordance with the product usage in this operation manual.
- Unplug the power plug from the outlet after use or before cleaning the product.
- The device has a heating function. Be careful when used by persons who are not sensitive to heat.
- Do not use on sluggish and overly sensitive skin or those with poor circulation part.
- Do not allow children or persons with physical and mental disabilities to use the warming function by themselves. Improper handling can be dangerous.
- To avoid fire, electric shock, or injury due to overheating, do not Operate the product under a pillow or blanket.
- Do not cover the use, so as not to cause the motor to overheat and cause malfunction or damage.
- The use of products by people with physical and mental disabilities and disabilities requires the guidance of others and can bear the responsibility responsible for its safety.
- If you feel uncomfortable during the massage, please stop the massager immediately.
- Do not use this product under the following conditions: The power cord or power plug is damaged, malfunctioning, dropped to the ground or damaged, dropped into water,Please contact the service center for inspection and repair.

- Do not use spare parts not recommended by the manufacturer.
- Do not insert or drop objects into any openings in the product.
- If you want to disconnect the power, turn off all operations, and then unplug the power plug from the outlet.
- This product is not intended for self-diagnosis or treatment. If diagnosis or treatment is required,It should be handled by professional medical personnel.

2.The use of the environment

- Do not use the product in a humid or dusty environment to avoid Fault.
- Do not use this product in an environment where the room temperature is higher than 40℃.
- Do not place the product outdoors.
- Do not place the product near electric heaters, near stoves, or in direct sunlight.


3.People who are not suitable to use this product

If you have the following conditions, please consult your doctor before use.


- Those who are being treated by a doctor or who are feeling unwell.
- Heart disease, cerebral arteriosclerosis, cerebral thrombosis.
- Patients with malignant tumor.
- Pregnant and menstruating women.
- Suffer from osteoporosis.

- Those using a pacemaker or implanted medical device.
  - Patients with acute illness.
  - Minors.
4. Pay attention
- Do not expose the control switch to water or liquids.
  - Do not operate the control switch with wet hands.
  - Do not fall asleep while using the product.
  - Do not use this product while intoxicated or unwell.
  - Do not use if the cloth cover of the product is damaged.
  - Do not use this product for more than 30 minutes continuously.
  - Please check the product for damage before use.
  - Please do not use a voltage other than the rated voltage.
  - Do not use one hour before or after eating to avoid adverse effects.

Product composition

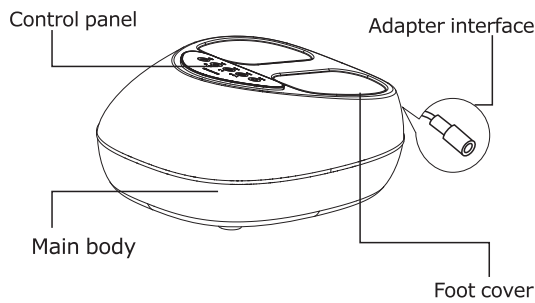


Removable foot cover

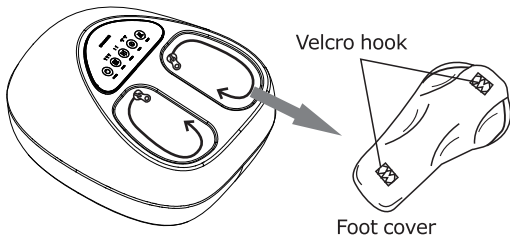

There is a zipper head here, which can be cleaned by opening the zipper and removing the foot cover. Please ensure that the foot cover is completely dry before putting it back into the foot massage machine.

Operating Instructions for Buttons

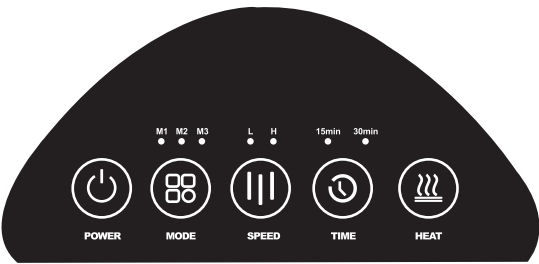

**POWER:** Click the  button to turn ON/OFF the foot massager  
In the standby state, click the power button, you can open the foot massager, at this time the Mode is M1, SPEED for the first gear, the hot compress function is open, the working time of 15 minutes.

**MODE:**Click the  button to switch airbag pressure mode  
M1: Low pressure airbag mode  
M2: Medium pressure airbag mode  
M3: High pressure airbag mode

**SPEED:** Click the  button to switch the massage speed  
L: Low Speed Massage  
H: High Speed Massage

**Time:**Tap the  button to switch the massage time  
15min: Automatically stops massaging after 15 minutes  
30min: Automatically stops massaging after 30 minutes

**HEAT:** Tap the  button to turn on/off the foot heat function  
When the word "HEAT" is seen, the heating function is on,if word is not displayed, the heating function is off.

Technical parameter

|               |                   |                |             |
|---------------|-------------------|----------------|-------------|
| Product name  | Foot massager     | Product model  | HXR-C060    |
| Input voltage | 110V-120V         | Rated voltage  | 12V         |
| Rated power   | 30W               | Product weight | About 2.6kg |
| Product size  | 375mm*332mm*178mm |                |             |

Troubleshooting

If you encounter problems when using this product,please refer to the following tips and possible solutions decision method. If this product is still in a faulty state, please send the product to a service center or Inspection and repair by authorized technicians.

| Abnormal situation       | Possible reason                                                                                                                                                                   | Approach                                                                                                                                                                                |
|--------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Not working              | 1. Not powered on.<br>2. The function keys are not turned on.                                                                                                                     | 1. Connect to the power supply or check whether the power supply is well connected.<br>2. Turn on the on/off button.                                                                    |
| Sudden stop of operation | 1. Power off.<br>2. The massage time exceeds 30 minutes.<br>3. Continuous use time is too long, Product overheating protection.<br>4. If the force is too large, press ball lock. | 1. Check whether the power interface is connected well.<br>2. Re-open the on/off button.<br>3. Re-open the on/off button after the product cools down.<br>4. Re-open the on/off button. |

Product keeper

1. Storage

- When not in use, please put the product in the color box and place it in a safe, dry, Cool and ventilated place.
- Do not twist the switch power supply and wires.
- Do not use switching power supplies and wires to hang the product.
- Avoid sharp objects scratching the surface of the product.
- Do not place the product in direct sunlight or high temperature.

2. Cleaning

- Before cleaning, make sure the main power switch is turned off and unplug it.
- Do not use switching power supplies and wires to hang the product.
- Never use aggressive cleaning agents, such as benzene or thinners to clean the product.
- Please use a dry cloth to clean the surface of the product.

3.Maintenance

- Keep the product clean and away from high humidity locations.
- If you do not use this product for a long time, please use a dust cover (please buy it yourself)Cover the product to prevent dust from falling into it.